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One-step
Emission
Technology

Gold Needle
Infusion
Technology


Qlin

ANION THERAPY COMB

PRODUCT MANUAL

address: 2-2, Incubator 3, Technology Park Malaysia, Bukit Jalil, 57000 Kuala Lumpur, Malaysia.

Introduction to Product Structure

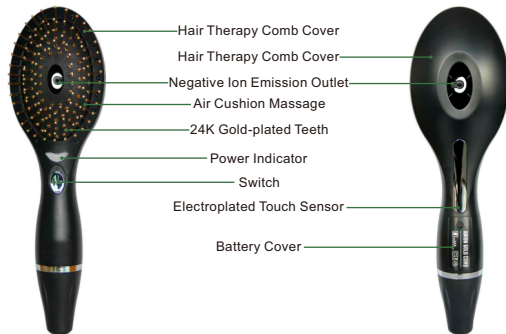Product Usage Instruction

This product must be inserted with two 7-cell batteries to achieve its unique function, please follow the steps below to install the battery properly:

1. Open the battery cover and insert two 7-cell batteries. Install the battery according to the positive and negative terminals, then replace the battery cover.
2. Slide the switch and the indicator lights up, showing that the super ion generator is starting to operate and actively releasing negative ions.
3. If you wish to use the maximum concentration of negative ions while combing your hair, make sure you grip the plated metal sensor on the handle upon using.
4. Upon completion of combing, slide the switch and notice that the indicator light is now off, a sign that the super negative ion generator has stopped working.

Special Note: When the switch is turned on and the indicator lights up, negative ions will be automatically generated. Any sound emitted from the negative ion generator is normal.

Special Features of the Product

1.This anion therapy comb is a hair comb purely derived by the principles of physics, resulting in a negative ion generator which boasts of the latest technology. It is completely free of any chemicals while it treats your hair, allowing you peace of mind while promising a 100% in perfect health care.

2.Subsequently, this anion therapy comb is suitable for both men and women of all ages, with its stylish appearance exclusively designed by American designers to add flair to a luxurious lifestyle. It has a built-in super ion emitter which has a life span of over 30,000 hours. The innovative negative ion technology allows this product to release a continuous flow of negative ions on its dual sides, granting full coverage of every inch of hair and scalp it combs through. It takes only one second to eliminate hair static, and another to smooth your dry hair, lock in moisture to repair damaged cuticles, so that your hair appears instantly silky smooth as you restore its lustre.

3.Negative ions are particles of negatively charged oxygen ions and water molecules in the air. They can penetrate into the hair core to suppress static electricity effectively while moisturizing hair, making hair easy to comb through as the negative ions penetrate into the hair layer. Every strand will turn supple as negative ions also work to combat germs, remove odour, relieve itch, lock in moisture, promote blood circulation, and etc. Eventually, after using the anion therapy comb for an extended period of time, your hair will become healthier naturally.

4.Simultaneously, negative ions can be helpful to our cerebral cortex. It can potentially assist the motor centre's function to regulate the subcortical sympathetic and parasympathetic balance of circulatory control.

5.In addition, negative ions can improve the ion transport mechanism of our vascular endothelial cells and restore the internal ion count within our cells. This is particularly important as based on research done by the Shanghai Institute of Blood Pressure, intracellular sodium increase is one of the important causes of hypertension and also a genetic trait of hypertension.

6.Furthermore, negative ions can unblock blood vessels by accelerating the "thrombosis recanalization". It can improve the condition of deformed red blood cells and prevent thrombosis within microcirculation. Moreover, it can help to prevent vasospasm and the detachment of hardened plaque fragments from the blood vessel walls.

Product Benefits

1.It is easy to carry, and it has a significant effect on easing toothache, migraine and otitis media.

2.The golden plated teeth will produce a high concentration of negative ions that can effectively remove aphids, activate facial cells, and produce beautifying and whitening results.

3.Negative ions of a small-diameter will be emitted in highly concentrated shots via this device. They will then penetrate through the blood-brain barrier to promote blood supply to the brain, strengthen the brain and calm nerves, aside from improving eye sight and providing quick relief to fatigue.

4.Neutralizes hair static, repairs scalp cells and damage to hair scales, making hair more supple, smooth and shiny.

5.Long-term use of this product can stimulate the acupuncture points of the head and enhance the regulating function of the cerebral cortex. It can also prevent hair loss effectively while promoting hair regeneration.

6.The negative ion emission outlet is close to the user's eyes. It can improve visual impairment such as myopia, glaucoma, and presbyopia, apart from combatting the viruses that cause refractory and trachoma.

7.The 24K gold-plated teeth of this comb stimulates body acupuncture points, effectively promoting blood circulation, improving health and the immune system. It is an essential tool for beauty, slimming and health therapy as it dredges the main and collateral channels of the body.

Double-sided Ion Release

360° Smoothing Function

Deep Locks
Moisture and Repairs

Seamless Air
Cushion Comb

Seamless Air
Cushion Comb

Acupoint Combing Method

Baihui Acupoint:

This point is a little way back from the top centre part of your head. You can comb from left to right, all the way down to this area.

Fengchi Acupoint:

It is 3 inches from the parallel point of your ear. You can comb from left to right, all the way down to this area.

Note: Combing through the entire head can help to repair damaged brain and hair cells. These two acupuncture points are often combed to prevent hair loss and strengthen hair roots.

Battery User Safety Precautions

1. Please use the No. 7 AAA dry carbon battery produced by regular manufacturers;
2. To extend the service life of this product, please use batteries of the same brand and same model.
3. Please do not use a mixture of batteries from different brands at the same time, such as old and new batteries, rechargeable and non-rechargeable batteries, or a combination of alkaline, carbon and zinc batteries.
4. Please make sure the batteries are installed according to the indicated positive and negative terminals.
5. Please replace the batteries immediately when the product indicator light shows that the batteries have weakened.
6. Do not remove the batteries from this product when it is not going to be in use for a long time.
7. Please dispose of your old batteries as per the relevant battery disposal regulations. All old batteries should be placed in the designated marked locations to prevent environmental pollution and reduce danger when recycling. Do not burn or bury your old batteries.

Warnings

The words of caution here are relevant to ensure personal safety and avoid property damage. Please make sure you adhere to the following.
Important: Do not immerse this device in water or other liquids.

- 1.This product should not be used, cleaned or maintained by children.
- 2.Please store this product in a place that is out of children's reach.
- 3.Do not use this product in wet areas or any place with water.
- 4.Please do not use this product together with any supplementary device that is not provided by the supplier.
- 5.Please inspect this product for damages prior to using it. If this product is already spoilt or shows signs of damage, please refrain from using it.
- 6.This product is not for the use of the disabled or those unable to use it under safe conditions without any supervision. Please use under guidance and strict surveillance to avoid mishaps.

Product Warranty Card

Your one year warranty starts from the day of your purchase

However, under the following conditions, there will be a small fee charged according to the standard price of all parts' or accessories.

- 1.Accidents caused by the carelessness of users (such as malfunctions or damage caused by movement or collision)
- 2.Product Failure or damage caused by force majeure (e.g. earthquake, tsunami, war and etc.)
- 3.Product failure and damage caused by the repair or modification of the product by any non-authorised after-sales service technician
- 4.Natural wear and tear of spare parts due to daily use
- 5.Product failure after the warranty period is over
- 6.Product damage due to improper packaging while taking it for repair
- 7.Product failure or damage caused by the transportation process (please contact the transportation department).

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